

# LEVO Kid & Junior

The every day wheelchair that let's  
your child stand, play and grow –  
everyday

LEVO Kid & Junior



LEVO stands for a better quality of life

**LEVO**.<sup>®</sup> **dt**   
THE EXPERTS IN STANDING

# No limitations for a life full of promise!



## To be part of...


Your child is active and wants to reach up and be part of it all. The LEVO Kid & Junior is their partner in play, in learning, in growing to their full potential. Like your child, the LEVO Kid & Junior is lively and fun with special abilities. Different seating positions and standing up to 86° offers a high degree of independence and flexibility. Easy adjustable growth plates, a wide range of accessories and options let the LEVO Kid & Junior smoothly grow with your child for many years.



# LEVO Stands for a better quality of life

## Highlights

- Smooth standing and positioning movement
- Compact, maneuverable, light
- Safe and secure stand-up operation elevating to any and all angles between sitting and 86° standing
- Switch activation and switch option for all levels of hand/finger function
- Adjustable for upper and lower leg length, foot position and back height
- Six ground stability points ensure solid and secure standing position
- User weight up to 265 lbs
- Wide range of different options and accessories
- Multiple color option
- Certificates: TÜV, CERAH, HI, FDA/510(k), Crash Test

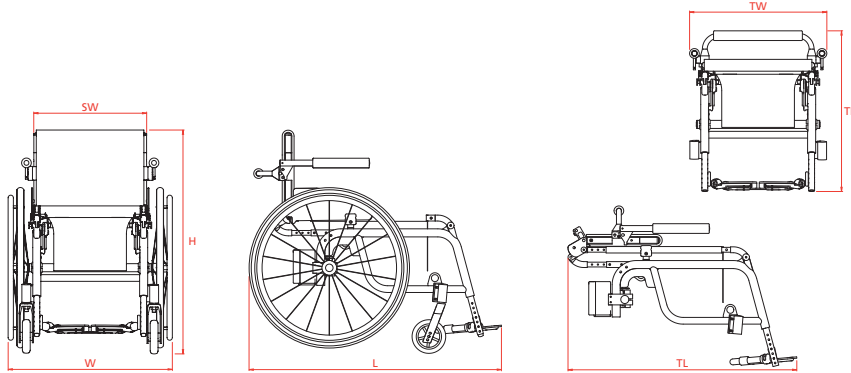

### Standard configuration

Armrests folding upwards and height and angle adjustable, folding down backrest angle adjustable, battery and charger, chest belt and knee support

### Options/accessories LCEV

Foam seat cushion, Anti-tip rollers, Tip aid, Backrest extension, Push handles, Push handle extension, Swing away chest support, Spinergy wheels, Camber, Spoke guard, V-Trak back, ... many more on request

## Technical Data

| Model                                | LCEV                             |           | Kid & Junior     |           |
|--------------------------------------|----------------------------------|-----------|------------------|-----------|
|                                      | S-M                              | L-XL      | Kid              | Junior    |
| Seat width (SW)                      | 13.8/15/15.7/16.5/17.7/18.9/19.7 |           | 11.8/13.8        |           |
| Total width (W)                      | 20.9/22/22.8/23.6/24.8/26/26.8   |           | 19.3/21.3        |           |
| Total length (L)                     | 33-35.4                          | 36.6-40.9 | 27.1-28.3        | 30-32.7   |
| Total height (incl. Backrest) (H)    | 31.5                             |           | 21.6             | 26        |
| Backrest upholstery height min.      | 12.2                             |           | 11               |           |
| Seat upholstery height front         | 19/20.1/21.2/22.4                |           | 13               | 15.7/17.3 |
| Seat upholstery depth                | 11.4-16.5                        | 17-22.4   | 9.4-15.7         |           |
| Footrest                             | Height and angle adjustable      |           |                  |           |
| Total weight (w/o detachable parts)  | 55 (42)                          | 57 (44)   | 37 (26)          | 42 (31)   |
| Rear wheels                          | 22"/24"/26"                      |           | 16"              | 20"/22"   |
| Batteries                            | 2 Batteries à 12V/2,7 Ah         |           |                  |           |
| Transport dimensions max. (TLxTWxTH) | 33x23.6x17.7                     |           | 24.4 x 17.7 x 17 |           |

Discrepancies in the measurements may arise due to wheelchair modifications. All measurements in inch and pound.

### Options/accessories Kid & Junior

Growth kit seat width 11.8 to 13.8 in, Growth kit Kid to Junior Height adjustable stroller handle, Attendant control for stand-up operation, Antitip devices, Table tray, Seat cushion, Back rest extension 2.4 to 6.3 in, V-Trak back, Spinergy wheels, ... many more on request

### Colors

5 standard colors (blue, red, yellow, platinum, black) and 190 other color options

